



Application Note AN-NIR-074

Surfactant in laundry detergent by Vis-NIR spectroscopy

Fast determination without using chemicals

Liquid laundry detergents contain fabric softeners, bleaching agents, surfactants, as well as enzymes. Out of these, the surfactant is the most important factor for the cleaning effect, as it breaks down the interface between polar and nonpolar compounds. This allows the detergent to be effective against greases as well as stains from soil or drinks.

Quantification of surfactant content is most

commonly performed by primary analyses (e.g., two-phase potentiometric titration). Disadvantages include manual sample preparation steps such as dilution and pH adjustment, and the method itself is time-consuming. In contrast, Vis-NIR spectroscopy has a time-to-result of less than 1 minute and does not require any sample preparation or chemicals for high quality data.

EXPERIMENTAL

A total of 37 samples with varying surfactant content were provided by a customer. The Vis-NIR spectra (Figure 2) were acquired on a Metrohm NIRS XDS RapidLiquid Analyzer equipped with 1 mm quartz cuvette (Figure 1). The samples were measured as-is, without any sample preparation steps. Data collection and model development was carried out with the Vision Air complete software package.



Figure 1. The NIRS XDS RapidLiquid Analyzer with a 1 mm quartz cuvette, used to collect the spectra of surfactant samples.

Table 1. Hardware and software equipment overview.

| Equipment | Metrohm number |
|--------------------------|----------------|
| XDS RapidLiquid Analyzer | 2.921.1410 |
| NIRS 1mm quartz cuvette | 6.7401.200 |
| Vision Air 2.0 Complete | 6.6072.208 |

RESULT

The obtained graph (Figure 3) displays a high correlation ($R^2 = 0.97$) between the values predicted by NIRS and the reference method. The nearly perfect

ratio of the SEC and SECV illustrates the validity of the model.

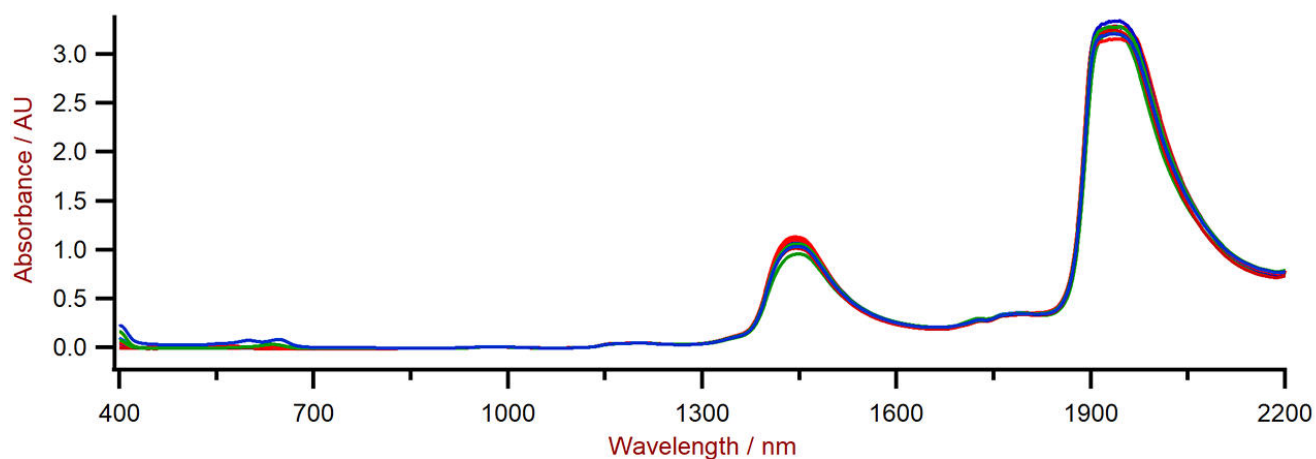


Figure 2 Selection of liquid detergent Vis-NIR spectra obtained using a XDS RapidLiquid Analyzer and a 1 mm quartz cuvette.

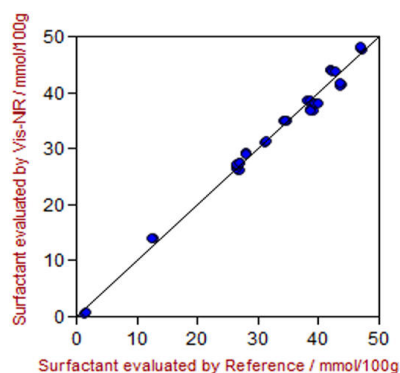


Figure 3 Correlation diagram and the respective figures of merit for the prediction of surfactant in liquid detergent using a XDS RapidLiquid Analyzer. The surfactant lab value was evaluated using HPLC.

Table 2. Figures of merit for the prediction of the surfactant content in liquid detergent using a XDS RapidLiquid Analyzer.

| Figures of merit | Value |
|------------------------------------|-----------------|
| R^2 | 0.97 |
| Standard error of calibration | 2.20 mmol/100 g |
| Standard error of cross-validation | 2.38 mmol/100 g |

CONCLUSION

The results presented herein show that the Vis-NIR method is excellently suited for the fast quantification of surfactant concentration in detergents. Using Vis-

NIR for this application saves 10 minutes per sample compared to other methods (Table 3).

Table 3. Time to result overview for the different parameters

| Parameter | Method | Time to result |
|----------------------|--------------------------|--|
| Surfactant (anionic) | Potentiometric titration | 10 min (adding solutions, stirring, pH-adjustments, determination) |

Internal reference: AW NIR CN-0015-102018

CONTACT

Metrohm USA
9250 Camden Field Pkwy
33578 Riverview, FL

info@metrohmusa.com



NIRS XDS RapidLiquid Analyzer

Rapid, precise analyses of liquids and suspensions of all types.

The NIRS XDS RapidLiquid Analyzer enables rapid, precise analyses of liquid formulations and substances. Precise measurement results at the push of a button make the NIRS XDS RapidLiquid Analyzer an equally reliable and simple solution for quality monitoring in laboratories and processes. The samples are transferred to quartz cuvettes designed for multiple use or disposable glass vials; a tempered sample compartment ensures reproducible analysis conditions and thus accurate measurement results.



Quartz cuvette, 1 mm

Quartz glass cuvette with a window made of quartz glass of maximum purity and homogeneity. The cuvettes have a transmission of more than 80% in the wavelength range of 200–2500 nm.

Compatible with:

- NIRS Spacer for set of 12.5 mm cuvettes (6.7403.180)
- DS2500 holder for 1 mm cuvette (6.7492.100)
- Holder OMNIS NIR, cuvette, 1 mm (6.07401.010)

Vision Air 2.0 Complete

Vision Air - Universal spectroscopy software.

Vision Air Complete is a modern and simple-to-operate software solution for use in a regulated environment.

Overview of the advantages of Vision Air:

- Individual software applications with adapted user interfaces ensure intuitive and simple operation
- Simple creation and maintenance of operating procedures
- SQL database for secure and simple data management

The Vision Air Complete version (66072208) includes all applications for quality assurance using Vis-NIR spectroscopy:

- Application for instrument and data management
- Application for method development
- Application for routine analysis

Additional Vision Air Complete solutions:

- 66072207 (Vision Air Network Complete)
- 66072209 (Vision Air Pharma Complete)
- 66072210 (Vision Air Pharma Network Complete)