



## Application Note AN-BAT-006

# Determination of the MacMullin number

## The stacking method

One important component of lithium ion batteries is the separator, a thin, electronically insulating foil which physically separates the two electrodes. Usually the separator is made of polymeric material.

In order to measure the performances of separators, the MacMullin number was introduced as the ratio of the ionic conductivity of the pure electrolyte, and the ionic conductivity of the separator foil filled with electrolyte.

In this application note, the MacMullin number is calculated by applying the so-called stacking method. This method consists of performing electrochemical impedance spectroscopy on samples with different stack thickness, achieved by increasing the number of separators inside the cell. Then the ion resistance is calculated from the data fitting and is plotted versus the number of separators in the cell. The slope gives the MacMullin number.

## INTRODUCTION

One important component of lithium ion batteries is the separator foil – a thin, electronically insulating, porous membrane that prevents the electrodes from being in direct contact with each other while being permeable to ionic current. To maintain a good lithium ion battery performance, the resistance of a separator foil filled with electrolyte should be low. Consequently, information about the effective ion conductivity of these systems is an important

parameter for the performance of a separator material.

Usually, the focus is not on the effective ionic conductivity, but on another quantity: the so-called MacMullin number  $N_M$  [1]. The MacMullin number is the ratio of the ionic conductivity of the pure electrolyte  $\sigma_{electrolyte}$  and ionic conductivity of the separator foil filled with electrolyte  $\sigma_{separator}$  (Equation 1).

$$N_M = \frac{\sigma_{electrolyte}}{\sigma_{separator}} \quad 1$$

Therefore, the MacMullin number is a quantity that describes the decrease of the effective conductivity by the presence of a separator foil. For lithium ion battery

separators, often a value between 4 and 20 is found [1–3]. Theoretically, the MacMullin number is linked to the separator foil's porosity  $\epsilon$  and tortuosity  $\tau$  [2, 3].

## EXPERIMENTAL

### Chemicals

As electrolyte, a 1 mol/L  $\text{LiPF}_6$  (lithium hexafluorophosphate) solution in a 1:1 ratio in volume mixture of EC (ethylene carbonate) and DMC (dimethyl carbonate, purchased from Sigma-Aldrich / Merck KGaA) was used without any further

purification. After receiving, the electrolyte was stored and handled inside of an argon-filled glove box. As battery separator material, a tri-layer separator based on polyethylene/polypropylene with a thickness of 21.5  $\mu\text{m}$  was used.

## Sample preparation and measuring setup

Circular specimens with a diameter of 12 mm were punched out from the separator foil. To ensure an optimal wetting, the specimens were stored for at least 24 hours in the electrolyte solution.


For the electrochemical measurements, a TSC battery standard cell (Figure 1) in combination with a Microcell HC setup was used.

As current collectors, two planar stainless steel disc electrodes (diameter 8 mm) press-fitted into a PEEK sleeve were used. The cell constant is the ratio of the electrode area  $A$  ( $\text{cm}^2$ ) and its thickness  $l$  (cm), and it is measured in  $\text{cm}^{-1}$ . To calculate the cell constant, the area of the stainless steel electrodes and the separator thickness were measured, resulting in a cell constant value of  $0.0043 \text{ cm}^{-1}$ .

The cell constant, together with the resistance related to the bulk ion transport inside of the porous separator network  $R_{ion}$  ( $\Omega$ ), is used in the calculation of the conductivity  $\sigma$  ( $\text{S cm}^{-1}$ , Equation 2).

$$\sigma = \frac{1}{R_{ion}} \cdot \frac{l}{A}$$

2

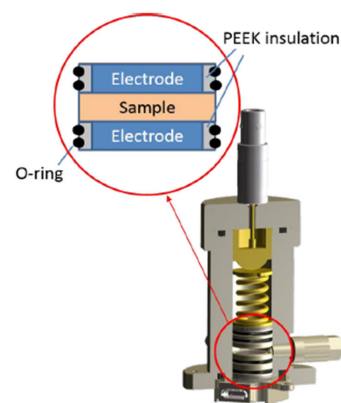

**Figure 1.** Schematic drawing of the measuring cell TSC battery standard. The separator foils filled with electrolyte were placed between the upper and lower electrode.

The contact pressure applied onto the electrolytesoaked separator foil stack inside of the TSC battery cell was adjusted to approximately 100 kPa using a gold-plated spring with a spring constant of  $2.43 \text{ N mm}^{-1}$ . At this contact pressure, no significant compression of the separator foil was found.

In the Microcell HC setup, the temperature is adjusted via a Peltier element, with an accessible temperature range of  $-40 \text{ }^\circ\text{C}$  to  $+100 \text{ }^\circ\text{C}$ . To accurately measure the temperature, a Pt100 temperature sensor is embedded in the base unit of the TSC battery cell, at a position very close to the sample. The accuracy of the

temperature measurement is  $0.1 \text{ }^\circ\text{C}$  with respect to the sensor position inside of the base unit. For the determination of the MacMullin number, the temperature was adjusted to  $20.0 \text{ }^\circ\text{C}$ .

Because the TSC battery measuring cell is completely airtight, the experiments took place on the Microcell system outside of the glove box.

For the impedance measurements, a Metrohm Autolab PGSTAT204 equipped with a FRA32M module was used. The communication with the temperature controller is integrated in the NOVA software enabling automated measuring routines.

## Measurement parameters

The impedance measurements were performed at a frequency range of 1 MHz to 100 Hz with an amplitude of 10 mV.

Impedance measurements were taken on samples with different stack thickness, achieved by increasing the number of separators inside the cell.

After reaching the temperature set point of 20.0 °C, a waiting time of 300 s was chosen to ensure for

complete thermal equilibrium.

**Table 1.** Steps and relative actions per step.


Step	Action to be performed
1	Transferring the cell into the glove box and loading the cell with one specimen of the electrolyte-filled separator foil.
2	Taking the prepared cell out of the glove box and connecting it to the Microcell HC cell stand.
3	Connecting the measuring device (2-electrode configuration).
4	Setting the temperature to 20.0 °C and waiting for 300 s.
5	Performing an impedance spectroscopy measurement.
6	After finishing step 5, repeating steps 1–5 by adding further specimens of the electrolyte-filled separator foil (one per loop) until a stack of five specimens is created.

## RESULTS

The resulting impedance spectra for different stack thicknesses are shown in **Figure 2**.

At high frequencies, the impedance behavior is dominated by inductive effects caused by the connection cables. The intersection at the  $Z'$ -axis is close to the bulk resistance value for bulk ion transport inside of the porous separator network. In absence of the inductive behavior at high frequencies, they would be identical. The increase at lower frequencies is then caused by electrode polarization.

As expected, increasing the stack thickness by adding further specimens of the electrolyte-soaked separator foil leads to a shift to higher values for the intersection point due to the higher bulk resistance value.

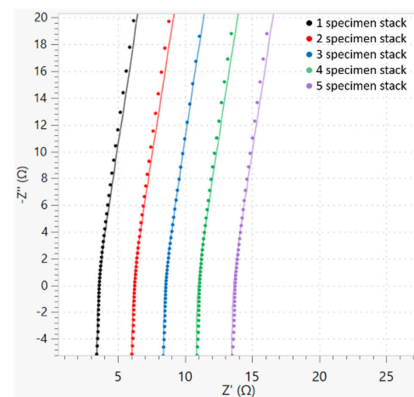

**Figure 2.** Impedance spectra of the soaked separator foils for different stack thicknesses. For a better visibility of the shift along the  $Z'$ -axis, the spectra were cut at 3 kHz. The points are the measured values and the lines represent the fit.

The spectra were fitted with the equivalent circuit shown in **Figure 3**.

The inductance  $l_{cable}$  takes into account the high frequency contributions by the electrical cables and connectors, the resistor  $R_{ion}$  represents the ion migration in the porous separator network, and the constant phase element  $CPE_{pol}$  describes the electrode polarization at low frequencies.

The resulting stack thickness-dependent values for  $R_{ion}$  are given in **Table 2**. The fit error of  $R_{ion}$  was lower than 0.4% in all cases.




**Figure 3.** Equivalent circuit for fitting the impedance spectra.

**Table 2.** Ion migration resistance per sample number.

$N_{specimen}$	$R_{ion}/\Omega$
1	3.38
2	5.98
3	8.34
4	10.78
5	13.43

A plot of these values shows a linear relationship, where the slope gives a value of  $2.49 \Omega$  for the ionic resistance change per added specimen  $\Delta R_{ion} / \Delta N_{specimen}$  (**Figure 4**).

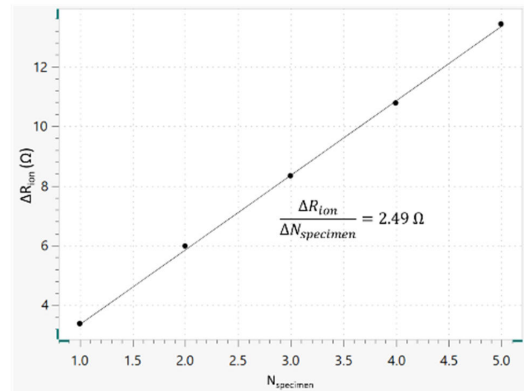

**Figure 4.** Plot of total ionic resistance value  $R_{ion}$  in dependence of the number of separator specimens  $N_{specimen}$ .

Taking into account the separator thickness  $d$ , the active area  $A$  of the stainless steel electrode, and the resistance per added specimen  $\Delta R_{ion} / \Delta N_{specimen}$ ,

the corresponding ionic conductivity of the electrolyte-soaked separator  $\sigma_{separator}$  can be calculated as shown in **Equation 3**.

$$\sigma_{separator} = \frac{1}{\Delta R_{ion}/\Delta N_{specimen}} \cdot \frac{d}{a} \quad 3$$

Here, a value of  $\sigma_{separator} = 1.7 \text{ mS cm}^{-1}$  was measured at 20 °C. The conductivity of the pure electrolyte solution at 20.0 °C was determined to be  $\sigma_{electrolyte} = 9.9 \text{ mS cm}^{-1}$ .

Therefore, the MacMullin number was calculated to be  $N_M = 5.8$  which lies within the typical range for such separators [3].

It should be mentioned that Raccichini et al. found that the MacMullin number can be slightly dependent on the chosen electrolyte solution [4], although

Landesfeind et al. showed that there was no dependence of the MacMullin number on the electrolyte salt, solvent composition, or salt concentration for the samples chosen by them [3]. Therefore, further research has to be performed to shed more light on the influence of the electrolyte properties on the determined MacMullin number values for separator foil. Ideally, the MacMullin number should be determined by using the electrolyte that is used for the final application.

## CONCLUSION

In this application note, the so-called stacking method is applied to determine the MacMullin number of separator foils. A Metrohm Autolab Microcell HC

setup combined with a TSC battery standard cell is used, together with a typical lithiumion battery electrolyte as well as a typical separator material.

## ACKNOWLEDGEMENT

The results presented here are part of a sub-project related to the project LiMES funded by the Federal Ministry for Economic Affairs and Energy as part of the

7. Energieforschungsprogramm which is highly appreciated.

## BIBLIOGRAPHY

1. P. Arora, Z. Zhang, Chem. Rev. 2004, 104, 4419-4462
2. M. J. Martínez, S. Shimpalee, J. W. Van Zee, J. Electrochem. Soc. 2009, 156 (1), B80-B85.
3. J. Landesfeind, J. Hattendorff, A. Ehrl, W. A. Wall, H. A. Gasteiger, J. Electrochem. Soc. 2016, 163 (7), A1373-A1387.
4. R. Raccichini, L. Furness, J. W. Dibden, J. R. Owen, N. García-Araez, J. Electrochem. Soc. 2018, 165 (11), A2741-A2749.

## CONTACT

Metrohm Suisse SA  
Industriestrasse 13  
4800 Zofingen

info@metrohm.ch

## CONFIGURATION



### Autolab PGSTAT204

Le PGSTAT204 associe faible encombrement et conception modulaire. Cet appareil comprend un potentiostat/galvanostat de base avec une tension disponible de 20 V et une intensité maximum de 400 mA ou 10 A en association avec le BOOSTER10A. Le potentiostat peut évoluer à tout moment au moyen d'un module complémentaire, comme le module de spectroscopie d'impédance électrochimique (SIE) FRA32M.

Le PGSTAT204 est un appareil d'un prix abordable qui trouve toujours une place dans le laboratoire. Il dispose d'entrées et de sorties analogiques et numériques pour contrôler les accessoires Autolab et les appareils externes. Le PGSTAT204 comprend un intégrateur analogique intégré. Associé au logiciel performant NOVA, il peut être utilisé pour la plupart des techniques d'électrochimie standard.



### Cellule de mesure TSC Battery Standard


La **TSC Battery** est une cellule de mesure fermée destinée à l'assemblage de **demi-cellules** ou de **cellules complètes de batterie** ou de **cellules à supercondensateur** avec, par ex., matières actives et séparateurs. Elle peut également être utilisée pour analyser des échantillons solides ou de type gel en contact avec des électrodes de métal dans une géométrie plan/plan.

Vous disposez d'une grande **flexibilité pour configurer la cellule**, c'est ainsi que vous pouvez effectuer **des mesures à deux et trois électrodes** par simple insertion sur le côté d'une micro-électrode de référence.

Le modèle standard comporte des électrodes de métal plates. Pour les collecteurs de courant conçus pour utiliser du lithium élémentaire comme matériau d'électrode, veuillez vous reporter au modèle étendu.

»

Cette cellule est hermétique, donc protégée contre l'air ambiant, et peut être utilisée avec des échantillons sensibles à l'humidité.

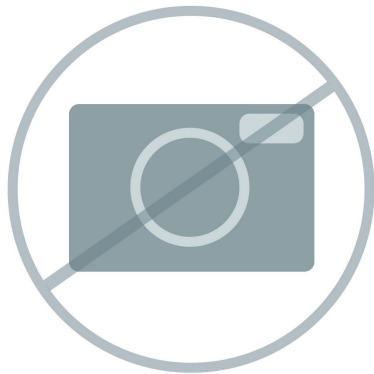

### Logiciel avancé pour la recherche électrochimique

NOVA est le logiciel conçu pour le contrôle de tous les instruments Autolab avec interface USB.

Conçu par des électrochimistes pour des électrochimistes, NOVA apporte plus de puissance et plus de flexibilité à votre potentiostat/galvanostat Autolab en intégrant plus de deux décennies d'expérience utilisateur et la toute dernière technologie logicielle .NET.

NOVA propose les fonctionnalités inédites suivantes :

- Un éditeur de procédures performant et flexible
- Une vue d'ensemble claire des données pertinentes en temps réel
- Des outils d'analyse de données et de tracés puissants
- Contrôle intégré des périphériques externes comme les instruments LQH Metrohm