



Application Note AN-RS-023

# Trace Detection of Aspartame in Beverages

## Protecting consumer safety with Misa

Aspartame is an artificial sweetener touted as a diet-conscious alternative to sugar in beverage and food products. Some studies suggest that consumption of aspartame is correlated with increased risk for brain and hematopoietic cancers, however, the majority of studies find aspartame to be a safe food additive. Consequently, the US and EU approve aspartame as a multi-purpose sweetener with an acceptable daily intake of 40 mg/kg body weight/day. However, the

clear health hazard to individuals suffering from phenylketonuria and ongoing criticism by health food advocates continues to fuel the challenge against aspartame's widespread use in the food industry. Using Misa (Metrohm Instant SERS Analyzer), beverage products are screened for aspartame levels with no sample preparation beyond simple dilution of a consumer product.

## INTRODUCTION


Misa is a versatile analytical tool for detecting additives in food products. This application note describes a facile method for the rapid and sensitive

detection of aspartame in carbonated water and diet cola.

## REFERENCE SPECTRUM AND LIBRARY CREATION

To establish a reference spectrum for aspartame, a pure standard dissolved in water at a concentration of 1 mg/g is analyzed using gold nanoparticles (Au NPs).

The unique baseline-corrected spectrum shown in **Figure 1** can be used to create a library entry for aspartame.

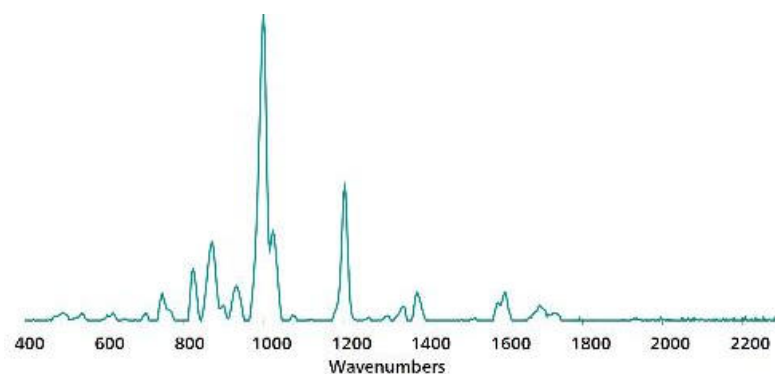

**Figure 1.** Standard SERS Au NP reference spectrum for aspartame in water.

## EXPERIMENT

A stock solution of aspartame in water was used to spike carbonated water from a local grocer to create a concentration range of samples: 100, 50, 10, 5, and 1 µg/mL, and 100 ng/mL. For testing purposes, 100 µL of each sample was pipetted into a glass vial containing 800 µL of Au NPs and 100 µL of 0.5 mol/L NaCl. The contents were shaken to mix, and the vial was placed into the Misa vial attachment for spectral acquisition. Diet cola from the same vendor was subjected to a 10x dilution (1:9 cola:water) and analyzed using the same procedure. Dilution reduces spectral contribution from other ingredients in the cola matrix.




**Table 1.** Experimental parameters

| Instrument           |                  | Acquisition |     |
|----------------------|------------------|-------------|-----|
| Firmware             | 0.9.33           | Laser Power | 5   |
| Software             | Misa Cal V1.0.15 | Int. Time   | 1 s |
| Misa Vial Attachment | 6.07505.040      | Averages    | 10  |
| ID Kit - Au NP       | 6.07506.440      | Raster      | ON  |

## RESULTS

Overlaid baseline-corrected SERS spectra acquired for carbonated water spiked with varying concentrations


of aspartame shows detection down to 100 ng/mL (Figure 2).

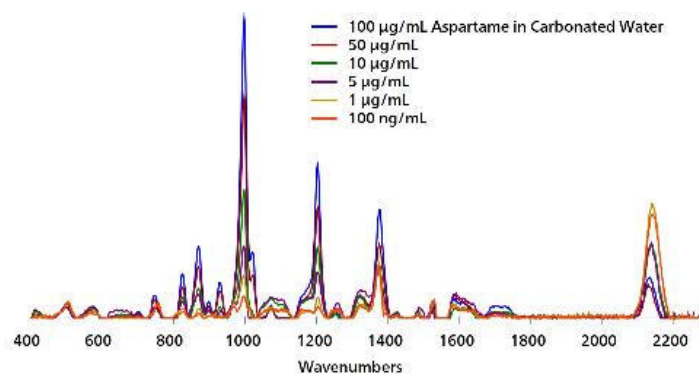

**Figure 2.** SERS Au NP concentration range for aspartame in carbonated water.

SERS analysis of 10x diluted diet cola with minimal processing provides high resolution detection estimated between 10–50 µg/mL, which is consistent

with typical concentrations of aspartame used to sweeten diet cola (e.g., 350 µg/mL).

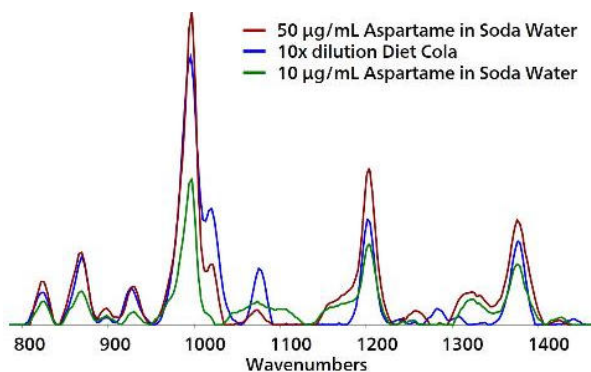

**Figure 3.** Detection of aspartame in 10x diluted diet cola.

## FIELD TEST PROTOCOL

### Detection of aspartame in the field

Using a pipette, add 1 drop of diet soda to a clean vial. Add 9 drops DI water to soda in vial, cap, and shake to mix. Fill a *clean vial* halfway full with Au NPs. Using pipettes, add 2 drops each of sample solution and

NaCl solution to Au NPs, then cap and shake the vial gently to mix. Insert into vial attachment on Misa for measurement.

**Table 2.** Requirements for field test protocol

|                |                            |
|----------------|----------------------------|
| ID Kit - Au NP | 6.07506.440                |
| includes:      | Gold nanoparticles (Au NP) |
|                | Scoop                      |
|                | Disposable pipettes        |
|                | 2 mL glass vials           |
| Reagents       |                            |
| DI water       |                            |
| NaCl solution  | 3 g NaCl in 100 mL water   |
| Test settings  | Use ID Kit OP on MISA      |

## CONCLUSION

The rapid and sensitive detection of aspartame in beverages is demonstrated using Misa and Au NPs. This assay demonstrates a quick and easy,

cost-effective and portable solution for detection of food additives in low-resource testing environments.

## CONTACT

Metrohm Nordic AS  
Sandviksveien 26  
1363 Høvik

[firmapost@metrohm.no](mailto:firmapost@metrohm.no)

## CONFIGURATION



### MISA Advanced

Metrohm Instant SERS Analyzer (MISA) is a high performance, portable analyzer system used for rapid, trace level detection / identification of illicit materials, food additives and food contaminants. MISA features a high-efficiency spectrograph equipped with Metrohm's unique Orbital-Raster-Scan (ORS) technology. It has a minimal footprint and extended battery life, perfect for on-site testing or mobile laboratory applications. MISA offers various Laser Class 1 attachments for flexible sampling options. Analyzer operation is available through Bluetooth or USB connectivity.

The MISA Advanced package is a complete package that allows the user to perform SERS analyses using Metrohm's nanoparticle solutions and P-SERS strips.

The MISA Advanced package includes a MISA Vial Attachment, a P-SERS Attachment, a ASTM Calibration Standard, a USB Mini Cable, a USB Power Supply and MISA Cal software for operating the MISA instrument. A ruggedized protective case is also provided to securely store the instrument and accessories.



### ID Kit – Au NP

The ID Kit - Au NP contains the components a Mira / Misa user requires to perform a SERS analysis using gold colloidal solution. The kit contains a disposable spatula, dropper, sample vials and a bottle of gold colloid.