

Application Note AN-PAN-1034

# Analysis of Bayer aluminate liquors using online thermometric titration

The Bayer Process is the method used to refine alumina from bauxite ore, as smelting aluminum directly from alumina is much more cost- and energy-effective. In this process, «aluminate liquors» are created by digesting the crushed bauxite with CaO and NaOH at high temperatures. Additionally, the CaO causticizes carbonate which forms in the alkaline solution from organic degradation and CO<sub>2</sub> absorption from the atmosphere. Contaminations are removed at various steps in the process, and the liquor is filtered from the alumina crystals before it is

recycled back to the digestion step. Before the spent liquor can be reused, a determination of the concentrations of the total hydroxide («caustic»), carbonate, and alumina is required.


This Process Application Note is focused on monitoring total hydroxide, carbonate, and alumina concentrations online in aluminate liquors via thermometric titration with either the 2060 TI Process Analyzer or the 2035 Process Analyzer - Thermometric from Metrohm Process Analytics.

## INTRODUCTION

Aluminum is used everywhere: in automobiles, bicycles, soft drink cans, cookware, and is even found in most antiperspirants, yet it does not occur in a natural state. Aluminum is a reactive base metal and is mainly refined from bauxite ore, which contains approximately 60% alumina ( $\text{Al}_2\text{O}_3$ ). To smelt aluminum directly from bauxite would be extremely costly due to its high melting point.

The Bayer Process was developed in the late 19th

century to extract alumina from bauxite, as purified alumina is much easier to smelt, and this cycle is still used by most alumina refineries today. The bauxite ore must be finely ground to increase surface area, and then mixed with cleaned spent liquor, lime ( $\text{CaO}$ ), and caustic soda ( $\text{NaOH}$ ). This slurry is digested at high temperatures under pressure for several hours (Figure 1).

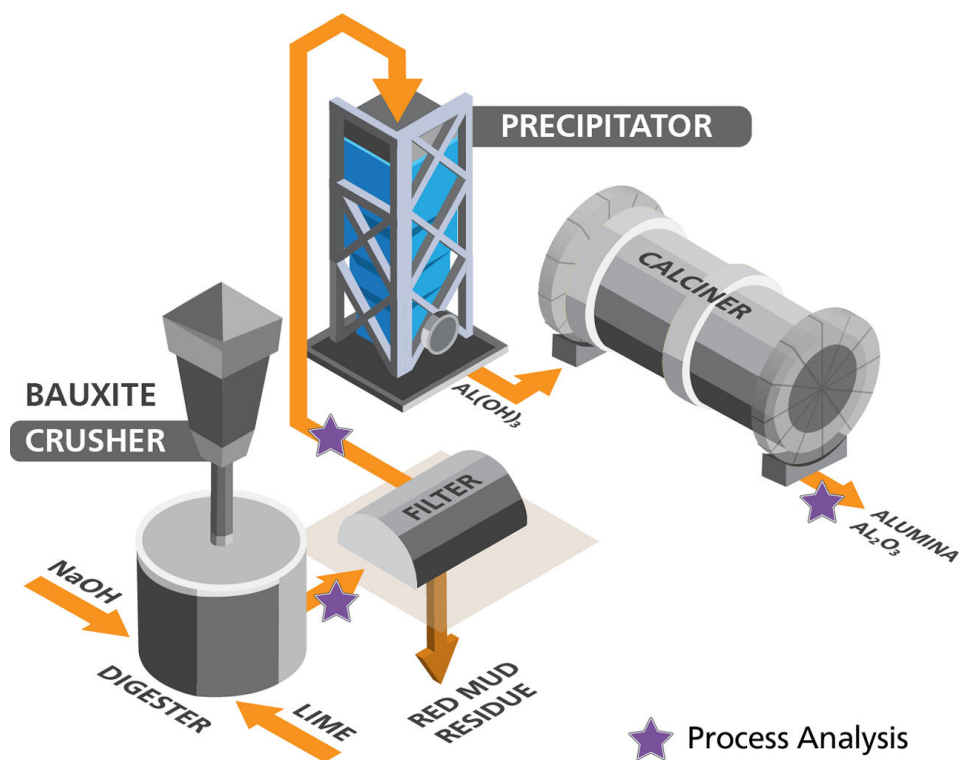

**Figure 1.** Bayer Process diagram with stars noting areas where online thermometric titration for process analysis can be integrated.




## APPLICATION

The sodium aluminate liquor is diluted with deionized water and complexed with sodium potassium tartrate, releasing one mole of hydroxide for each


The total hydroxide content of the liquor (total caustic) and the carbonate (total soda) content are

Potassium fluoride solution is then added to destroy the aluminotartrate complex, forming insoluble potassium sodium aluminum fluoride and releasing three moles of hydroxide (also determined by HCl) for each mole of aluminate (**Reaction 3**).

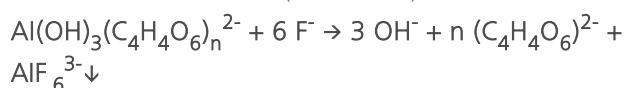


Reaction 3.

The sodium aluminate liquor is diluted with deionized water and complexed with sodium potassium tartrate, releasing one mole of hydroxide for each mole of aluminate present (**Reaction 1**).

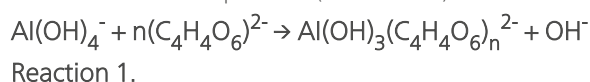


The total hydroxide content of the liquor (total caustic) and the carbonate (total soda) content are determined by titration with HCl (**Reaction 2**).

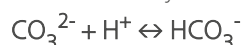

Reaction 2.

Potassium fluoride solution is then added to destroy the aluminotartrate complex, forming insoluble potassium sodium aluminum fluoride and releasing

three moles of hydroxide (also determined by HCl) for each mole of aluminate (**Reaction 3**).

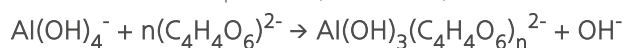


Reaction 1.

determined by titration with HCl (**Reaction 2**).  $\text{CO}_3^{2-} + \text{H}^+ \leftrightarrow \text{HCO}_3^-$  Reaction 2.

A second titration is then automatically and immediately performed to determine the aluminate content (as «alumina»). Total caustic is defined as the total hydroxide content of the liquor comprising unassociated hydroxide ions, and one hydroxide ion of the four found in the aluminate  $[\text{Al(OH)}_4^-]$  anion. Total soda is defined as the sum of the total caustic content plus the carbonate content of the liquor.

three moles of hydroxide (also determined by HCl) for each mole of aluminate (**Reaction 3**).

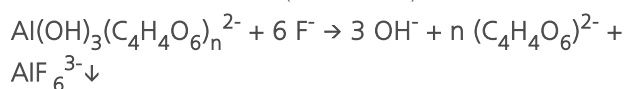

Reaction 3.

A second titration is then automatically and immediately performed to determine the aluminate content (as «alumina»). Total caustic is defined as the total hydroxide content of the liquor comprising unassociated hydroxide ions, and one hydroxide ion of the four found in the aluminate  $[\text{Al(OH)}_4^-]$  anion. Total soda is defined as the sum of the total caustic content plus the carbonate content of the liquor.





**Figure 3.** 2060 TI Process Analyzer from Metrohm Process Analytics.



**Figure 4.** 2035 Potentiometric Analyzer - Thermometric.

**Table 1.** Different parameters measured online with thermometric titration during Bayer Process monitoring.

Parameters	Range
Total caustic	17–150 g/L (as Na <sub>2</sub> O)
Total soda	1–155 g/L (as Na <sub>2</sub> O)
Alumina	17–170 g/L (as Al <sub>2</sub> O <sub>3</sub> )

## REMARKS

Highly concentrated liquors may need a reduced sample size and modified titrant quantities to effectively complex all aluminate with the tartrate reagent. Very dilute liquors may be titrated directly.

Pure sodium aluminate solutions are also produced for use in water purification, the manufacture of paper and of synthetic zeolites; the method described here is also suitable for these solutions.

## CONCLUSION

The 2060 TI Process Analyzer and 2035 Process Analyzer - Thermometric from Metrohm Process Analytics can not only measure the concentration of alumina, but also the total hydroxide and carbonate concentration in aluminate liquors via thermometric

titration. This method is the preferred solution since it is suitable for aggressive matrices, does not require sensor maintenance, and is a highly sensitive analysis technique.

## RELATED APPLICATION NOTES

[AN-PAN-1037 Online measurement of the acid number \(AN\) in oils with thermometric titration](#)  
[Brochure: 2060 TI Process Analyzer – Maximum flexibility for the toughest challenges in process](#)

[analysis](#)  
[Brochure: 2035 Process Analyzer – Multi-purpose analyzer for the online monitoring of industrial processes and waste waters](#)

## BENEFITS FOR THERMOMETRIC TITRATION IN PROCESS

- Detect process upsets via automated analysis
- Increased product throughput, reproducibility, production rates, dosage of chemicals, and profitability
- Fully automated diagnostics – automatic alarms for when samples are out of specification parameters



## CONTACT

Metrohm Portugal  
R. Frei Luis de Granada 14G  
1500-680 Lisboa

[vendas@metrohm.pt](mailto:vendas@metrohm.pt)

## CONFIGURATION



### 2035 Process Analyzer - Thermometric

The 2035 Process Analyzer for Thermometric Titration performs automated titrations with a fast-responding, highly sensitive temperature sensor. Instead of the electrochemical potential, the endpoint is determined by enthalpy, i.e., by change in temperature of the solution during the titration. Thermometric titration is by far the most robust titration method available on the market and is excellent for 24/7 online process applications, such as the monitoring of etching baths. No sensor calibration is necessary, and less cleaning steps are needed. Rapid analysis is possible with this technique; acid mixtures, e.g., can be analyzed in less than three minutes.

Thermometric titration can be used for a wide variety of titration analyses and is well-suited to handle aggressive sample matrices because of the robust thermometric sensor. The sensor needs no maintenance because fouling and other undesired interactions are minimal, and there are no membrane or diaphragm problems as with other titration methods. Thermometric titration is a problem solver for difficult samples which cannot be titrated potentiometrically, and is also a preferred technique in situations when HF is present in samples.



## 2060 Process Analyzer

The 2060 Process Analyzer is an online wet chemistry analyzer that is suitable for countless applications. This process analyzer offers a new modularity concept consisting of a central platform, which is called a «basic cabinet».

The basic cabinet consists of two parts. The upper part contains a touch screen and an industrial PC. The lower part contains the flexible wet part where the hardware for the actual analysis is housed. If the basic wet part capacity is not sufficient enough to solve an analytical challenge, then the basic cabinet can be expanded to up to four additional wet part cabinets to ensure enough space to solve even the most challenging applications. The additional cabinets can be configured in such a way that each wet part cabinet can be combined with a reagent cabinet with integrated (non-contact) level detection to increase analyzer uptime.

The 2060 process analyzer offers different wet chem techniques: titration, Karl Fischer titration, photometry, direct measurement and standard additions methods.

To meet all project requirements (or to meet all your needs) sample preconditioning systems can be provided to guarantee a robust analytical solution. We can provide any sample preconditioning system, such as cooling or heating, pressure reduction and degassing, filtration, and many more.