

# Identification of Forensic Fabrics Using a Portable Raman Spectrometer

At a crime scene, a police officer collects a fiber sample that may prove to be invaluable evidence in identifying a criminal or exonerating an innocent person. While FTIR has been used for analysis in the past, the strong absorption of the fabric or the glass slide where it is mounted makes the spectrum very hard to interpret. In recent years, Raman spectroscopy

has been studied extensively for forensic fiber analysis because of the high selectivity of Raman signatures, non-destruction nature of the test, and the ability to conduct the analysis without any sample preparation. The Raman test can be performed directly on fabrics or fibers mounted on glass slide with very little interference from the mounting resin or the glass.

## EXPERIMENTAL

In this study, six types of undyed fabric samples were tested: diacetate, bleached cotton, polyester, polyamide (nylon), acrylic, and wool. A B&W Tek i-Raman EX portable Raman spectrometer with 1064 nm laser excitation along with a fiber optic probe holder was used. The video microscope sampling accessory can be used for testing on thin fibers and is very useful for microsamples, and for looking at specific spots on a sample.

The identification test involves creating a library and comparing each individual fabric spectrum with the spectra in the library. The software BWID was used to generate the library as well as to conduct the identification. The resulting "Match" or "No Match" is based on HQI (hit quality index), which measures the level of correlation of the sample spectrum against a reference spectrum calculated by using the correlation coefficient algorithm. The HQI threshold for a "Match" result was set at 80, which indicates an 80% correlation score between the sample spectrum and the reference spectrum.

Four out of the six types of fabrics are unambiguously identified. The overlaid spectra for these four fabrics are shown in **Figure 1**. As demonstrated in the diagonal line in the comparison results of **Table 1**, there is clear differentiation to separate diacetate, bleached cotton, polyamide, and acrylic fabrics.



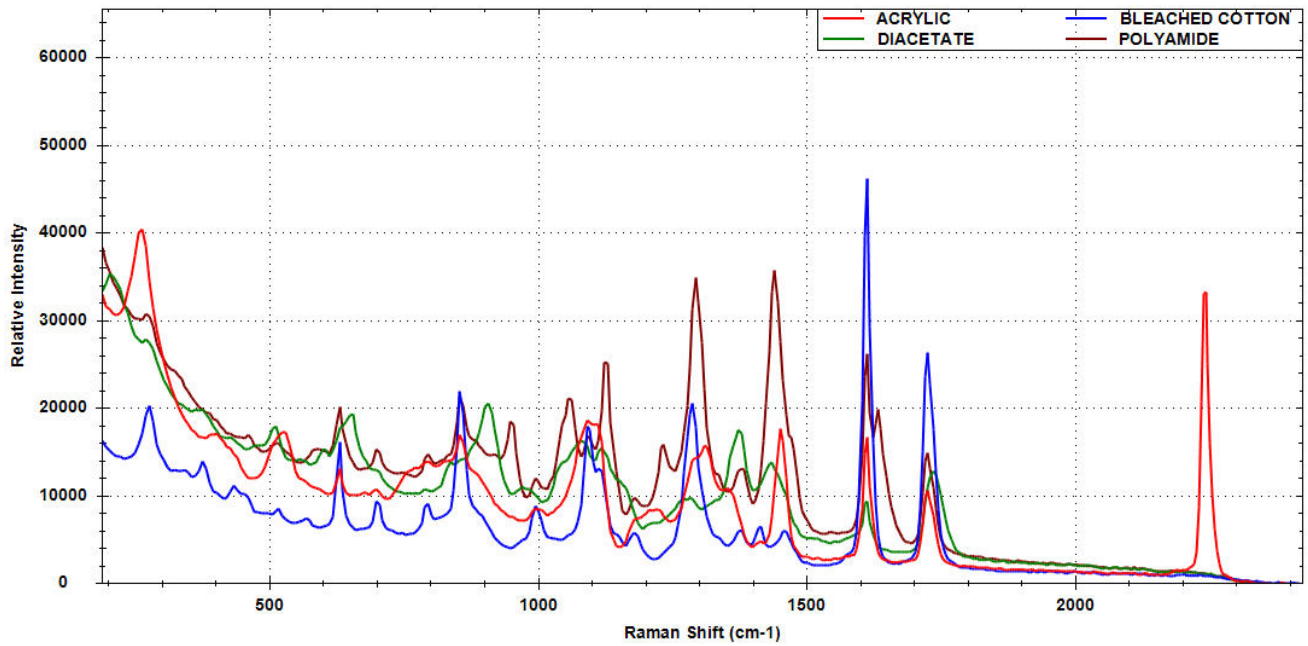


Figure 1. Overlay of spectra for acrylic, bleached cotton, polyamide, and diacetate

Library Sample \	Diacetate	Bleached Cotton	Polyamide	Acrylic	Polyester	Wool
Diacetate	Match (HQI=87.68)	No match	No match	No match	No match	No match
Bleached Cotton	No match	Match (HQI=94.08)	No match	No match	No match	No match
Polyamide	No match	No match	Match (HQI=91.98)	No match	No match	No match
Acrylic	No match	No match	No match	Match (HQI=93.97)	No match	No match
Polyester	No match	No match	No match	No match	Match (HQI=96.59)	2 <sup>nd</sup> Match (HQI=81.23)
Wool	No match	No match	No match	No match	2 <sup>nd</sup> Match (HQI=94.75)	Match (HQI=85.38)

Table 1. Identification results using BWID software

Polyester and wool are difficult to differentiate by HQI, as the Raman spectra are highly similar (Figure 2). However, since fibers made from animal hair contain protein keratin, an amide I band from the amino acid cysteine in the region from 1600-1690  $\text{cm}^{-1}$  would be expected[1] for wool and not polyester. In addition, since cysteine provides the disulphide crosslinks that

hold the polymer chains in wool, a disulphide S-S band at 523  $\text{cm}^{-1}$  would also be expected[1]. These two peaks that are distinctively related to animal proteins can be seen in the wool spectrum shown in Figure 2, with the amide I band at 1653  $\text{cm}^{-1}$  and the S-S band at 523  $\text{cm}^{-1}$ . These two unique peaks can be used to differentiate wool from polyester.

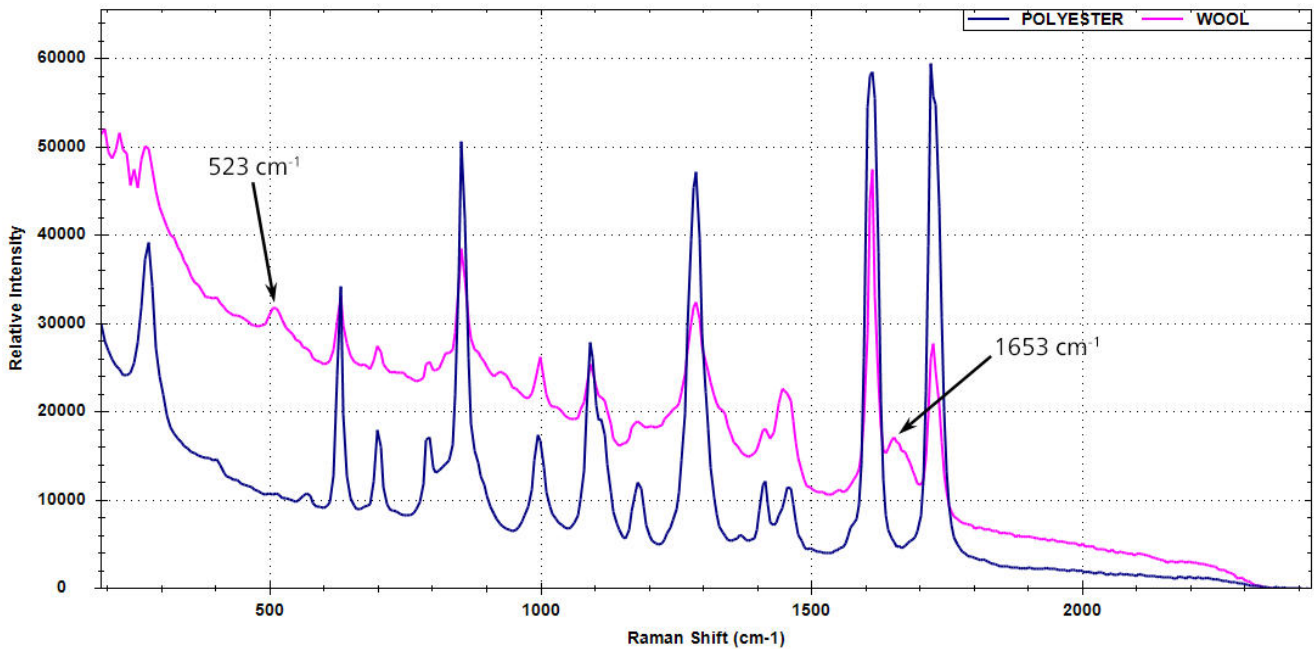


Figure 2. Overlay of spectra for wool and polyester

## CONCLUSION

With unique discrimination power, Raman spectroscopy is a powerful technology that can be applied to forensic fabric and fiber analysis.

Identification of an unknown fabric is achieved in several minutes, making it a practical choice for rapid identification either on site or in the forensic lab.

## REFERENCE

1. Li-Ling Cho. Identification of textile fiber by Raman microspectroscopy. J Forensic Science 2007; 6 (1):55-62

## CONTACT

Metrohm Viet Nam  
Phan Dinh Giot  
70000 Ho Chi Minh

info@metrohm.vn


## CONFIGURATION



### i-Raman EX Portable Raman Spectrometer

The i-Raman<sup>®</sup> EX is part of our award-winning series of i-Raman portable Raman spectrometers with our patented CleanLaze<sup>®</sup> laser with 1,064 nm laser excitation. Using a high-sensitivity InGaAs array detector with TE deep cooling, high dynamic range, and a high throughput spectrograph design, this portable Raman spectrometer delivers a high signal-to-noise ratio without inducing autofluorescence, making it possible to measure a wide range of natural products, biological samples (such as cell cultures), and colored samples.

The i-Raman EX provides a spectral coverage range from 100 cm<sup>-1</sup> to 2,500 cm<sup>-1</sup>, enabling you to measure across the entire fingerprint region. The system's small footprint, lightweight design, and low power consumption ensure research-grade Raman analysis capabilities at any location. The i-Raman EX comes equipped with a fiber probe and an XYZ-positioning-stage with probe holder. It can be used with a range of sampling accessories to facilitate measurements on a wide range of different samples. For expanded analysis capabilities, it can be used with our proprietary BWIQ<sup>®</sup> multivariate analysis software or BWID<sup>®</sup> identification software. With the i-Raman EX, you always have a high precision Raman solution for qualitative and quantitative analysis without fluorescence.

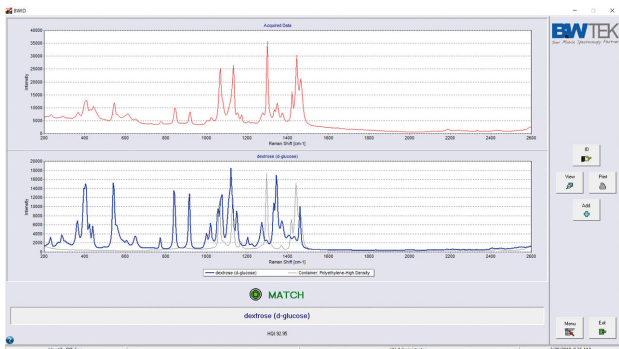

### BWID Standard Raman ID Software

Data acquisition, material identification and reporting software for B&W Tek portable Raman spectrometers to rapidly identify or verify materials based on user-built libraries, proprietary B&W Tek and third-party libraries. Includes system performance tests. One user license provided. Free upgrades included within one year of purchase date.



### Raman Video Micro-Sampling System (1064 nm)

Video microscope sampling system for use with B&W Tek's lab and industrial Raman probes. Includes a 20x objective at a working distance of 16 mm. Offers manual rough and fine adjustment on the X, Y, and Z-axes, coaxial LED illuminator for target alignment, video camera for sample observation, and is compatible with standard microscope objectives. Probe not included, available separately. 1,064 nm configuration.

BAC151C-1064