

## 1 Purpose

This document describes the improvements in the following new firmware versions in comparison with the previous versions:

- 5.779.0012 (replaces 5.779.0011)
- 5.779.0022 (replaces 5.779.0021)

## 2 Improvements

Topic	Description
Automatic recognition	The Windows Device Manager recognizes Optrodes as soon as the USB plug is plugged in.
Changing the wavelength	Starting with tiamo 2.5, the wavelength can be changed in the <b>DET U</b> , <b>MET U</b> and <b>MEAS U</b> commands. It is not possible anymore to change the wavelength with a magnet during measurements.


## 3 Compliance

The current firmware versions do not contain any modifications that affect the conformity of the instruments regarding 21 CFR Part 11.

**Herisau, September 15, 2015**



P. Hunziker  
Vice President,  
Head of Development

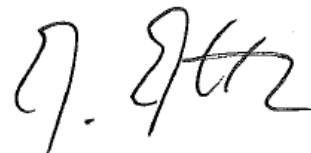

M. Matter  
Head of Quality Management